


# How to install the Pathomation Fiji plugin

Fiji plugin installation (Figure 1: how to install the plugin):

1. Extract the "HistoJ\_2.0.0.02c41a565e\_RC1.zip" file. You can get the plugin [here](#). Open the Fiji application and expand the "Plugins" menu by clicking on it.
2. Click on the "Install..." menu item.
3. Browse the extracted "HistoJ\_2.0.0.02c41a565e\_RC1.jar" file and click the "Open" button.
4. This will prompt another window to save the plugin. Confirm the installation by clicking the "Save" button.
5. Please click the "OK" button and restart the Fiji application to complete the installation.

Figure 1: how to install the plugin

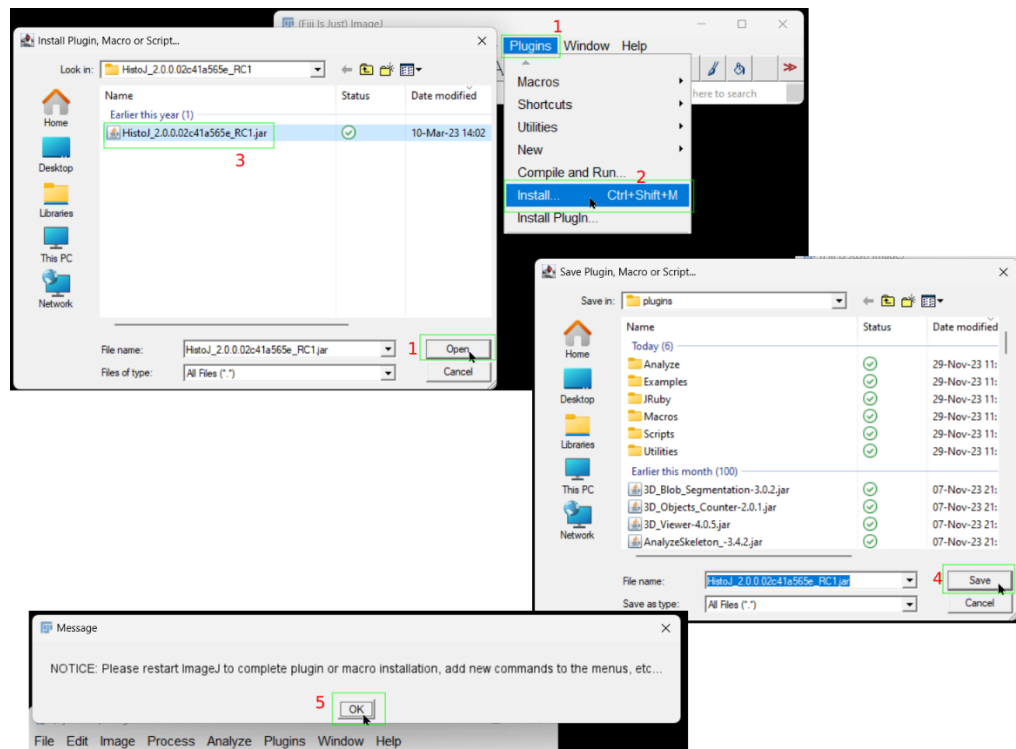

Fiji plugin startup guide (Figure 2: how to start):

6. Run PMA.start. You can get the latest version of PMA.start [here](#). Expand the "Plugins" menu by clicking on it.
7. Click the "HistoJ" plugin.
8. Click the "PMA.start" button.
9. Browse the image you want to use in Fiji.
10. Confirm your choice by clicking on the "OK" button.
11. Zoom into the preferred ROI. You can also leave the image as is. Using a "Small" scale in the next step will diminish the zoomed-in quality. Confirm your choice by clicking the "OK" button.
12. In the "Choose scale" window, you can choose the "Scale" of the image. Note that the larger you choose the scale to be, the longer the Fiji application will take to render the image.
13. Confirm your choice by clicking the "OK" button.

Figure 2: how to start

